


RESEARCH ARTICLE

OCT-angiography: A qualitative and quantitative comparison of 4 OCT-A devices

Marion R. Munk^{1,2*}, Helena Giannakaki-Zimmermann^{1,2}, Lieselotte Berger^{1,3}, Wolfgang Huf⁴, Andreas Ebner^{1,3}, Sebastian Wolf^{1,2,3}, Martin S. Zinkernagel^{1,2,3}

1 Department of Ophthalmology, Inselspital, Bern University Hospital, University of Bern, Switzerland, **2** Bern Photographic Reading Center, University of Bern, Switzerland, **3** Department of Clinical Research, Inselspital Bern University Hospital, University of Bern, Switzerland, **4** Center for Medical Physics and Biomedical Engineering, Medical University of Vienna, Vienna, Austria

* marion_munk@hotmail.com

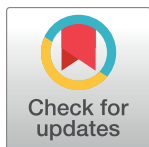

Abstract

Purpose

To compare the quality of four OCT-angiography(OCT-A) modules.

Method

The retina of nineteen healthy volunteers were scanned with four OCT-devices (Topcon DRI-OCT Triton Swept-source OCT, Optovue RTVue-XR, a prototype Spectralis OCT2, Heidelberg-Engineering and Zeiss Cirrus 5000-HD-OCT). The device-software generated en-face OCT-A images of the superficial (SCP) and deep capillary plexuses (DCP) were evaluated and scored by 3 independent retinal imaging experts. The SCP vessel density was assessed using Angiotool-software. After the inter-grader reliability assessment, a consensus grading was performed and the modules were ranked based on their scoring.

Results

There was no significant difference in the vessel density among the modules (Zeiss $48.7 \pm 4\%$, Optovue $47.9 \pm 3\%$, Topcon $48.3 \pm 2\%$, Heidelberg $46.5 \pm 4\%$, $p = 0.2$). The numbers of discernible vessel-bifurcations differed significantly on each module (Zeiss 2 ± 0.9 bifurcations, Optovue 2.5 ± 1.2 , Topcon 1.3 ± 0.7 and Heidelberg 0.5 ± 0.6 , $p \leq 0.001$). The ranking of each module differed depending on the evaluated parameter. In the overall ranking, the Zeiss module was superior and in 90% better than the median (Bonferroni corrected p -value = 0.04). Optovue was better than the median in 60%, Topcon in 40% and Heidelberg module in 10%, however these differences were not statistically significant.

Conclusion

Each of the four evaluated OCT-A modules had particular strengths, which differentiated it from their competitors.

OPEN ACCESS

Citation: Munk MR, Giannakaki-Zimmermann H, Berger L, Huf W, Ebner A, Wolf S, et al. (2017) OCT-angiography: A qualitative and quantitative comparison of 4 OCT-A devices. PLoS ONE 12(5): e0177059. <https://doi.org/10.1371/journal.pone.0177059>

Editor: Knut Stieger, Justus Liebig Universitat Giessen, GERMANY

Received: March 14, 2017

Accepted: April 23, 2017

Published: May 10, 2017

Copyright: © 2017 Munk et al. This is an open access article distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited.

Data Availability Statement: All relevant data are within the paper and its Supporting Information files.

Funding: The authors received no specific funding for this work.

Competing interests: I have read the journal's policy and the authors of this manuscript have the following competing interests: Sebastian Wolf is a consultant for Heidelberg Engineering and Zeiss. Martin S Zinkernagel is a consultant for Heidelberg Engineering. The remaining authors have no

financial/proprietary interests relevant to this paper. This does not alter our adherence to PLOS ONE policies on sharing data and materials.

Introduction

OCT-angiography (OCT-A) is a new diagnostic tool and heavily promoted as an alternative or an adjunct to classic fluorescein angiography (FA). It is a fast imaging tool, detecting streaming blood, thereby allowing to construct an image of the retinal vasculature; in contrast to “classical” FA it is dye-free, and therefore lacks significant side effects associated with the fluorescein injections such as such as vomiting, hypersensitivity reactions and cardiovascular complications. [1] This new technology allows the in situ, high-resolution visualization of the individual vascular layers. In contrast to FA, which displays only the superficial capillary network, OCT-A visualizes the superficial, the deep and the choroidal vascular network; even the middle capillary plexus can be identified. [2].

Several OCT manufacturers now offer OCT devices including algorithms enabling the practitioner to obtain regular OCT B-scans as well as volumetric angiographic images. Different techniques such as Doppler shift, speckle variance/decorrelation, phase variance, optical micro-angiography and correlation mapping are employed to differentiate blood vessels by depicting the change in the OCT-signal induced by the moving blood cells. [3,4] So far Angiovue optical coherence tomography angiography (Optovue RTVue XR Avanti, Optovue Inc., Fremont, CA) based on a split spectrum amplitude decorrelation angiography algorithm (SSADA), Zeiss AngioPlex (Cirrus HD-OCT 5000, Zeiss Meditec. Inc.) based on micro-angiography (OMAG) and SS-OCT Angiography employed in a Swept source OCT DRI OCT Triton (Topcon DRI OCT Triton Swept source OCT, Topcon, Japan) using the so called OCT angiography ratio analyses (OCTARA) algorithm are commercially available. Prototypes of the Spectralis OCT2 module (Heidelberg Engineering, Germany) with a full spectrum amplitude decorrelation algorithm, and a prototype of AngioScan (RS-3000 Advance OCT, Nidek Co., Ltd., Japan) based on a complex decorrelation algorithm were introduced and are currently tested. Further there are other OCT-A modules under development such as the OCT-A system inbuilt in the Copernicus Revo and REVO NX by OPTOPOL.

Previous studies aimed to compare the performance of different OCT-A techniques applied in the above listed OCT-A modules including phase variance, absolute complex difference, speckle variance and absolute intensity difference. It was confirmed that all methods generate excellent flow motion contrast images, with phase variance and absolute complex difference methods requiring more complex analyses than intensity based algorithms such as speckle variance and absolute intensity difference. [5] A recent study by De Vitis et al. compared the AngioVue (Optovue) with the Angioplex (Zeiss) and found that the Angioplex required shorter acquisition time and showed a lower number of motion artefacts when compared to the Angiovue. [De Vitis, 2016 #881].

However so far no data are available which systematically compare the commercially available OCT-A modules.

The aim of this study therefore was to compare the quality of 4 different commercially available OCT-A modules in healthy subjects.

Methods

Patients and setting

Nineteen healthy subjects were evaluated in this cross sectional study. Subjects had a visual acuity of 20/20 or better without a clinical history and without any evidence of an eye disease including retinal disease or glaucoma. Exclusion criteria were also the presence of diabetes or hypertension or any other cardiovascular disease. The retina was scanned using a Zeiss Cirrus 5000 HD-OCT (Zeiss Meditec. Inc, Germany), an Angiovue, RTVue XR Avanti (Optovue,

Inc), a Topcon DRI OCT Triton Swept source OCT and a prototype of Spectralis OCT2 (Heidelberg Engineering, Germany). The Topcon DRI Swept source (SS)- OCT used a wavelength of 1050nm, whereas the remaining devices used shorter wavelengths around 800nm (870nm, Heidelberg; 840nm, Optovue and Zeiss). The Zeiss Cirrus HD-OCT Model 5000 with Angio-plex uses a so called OCT- microangiography complex algorithm (OMAG) and an A- scan rate of 68Khz. OMAG identifies changes in the phase and intensity information of the OCT scans to quantify motion contrast.[6] For eye tracking the FastTrac technology is implemented and the retina is sampled 15 frames per second to minimize motion artefacts. Only areas which may be affected by motion artefacts are rescanned, which decreases the acquisition time. A 3x3 pattern with a 245x245 resolution was chosen, with a mean distance of 12.2 microns between each scan and each B-scan was repeated 4 times at the same position. The A scan depth is 2mm with an axial resolution of 5 μ m and a transverse resolution of 15 μ m.[6] The Optovue Angiovue utilizes SSADA, which splits the spectrum into different, small bands while employing a decorrelation measure. [7] With Optovue, a 3x 3-volume scan centered on the fovea was obtained with an A-scan rate of 70kHz. Each volume scan consists of 304x304 A-scans with 2 consecutive B-scans at each position. Two right angled OCT-A volumes scans are performed for orthogonal registration to correct for motion artefacts. [7] OCTA Ratio Analyses (OCTARA) employed by Topcon is an intensity ratio analyses and is not based on amplitude decorrelation. It does not require splitting the spectrum and therefore preserve axial resolution, which is important as SS-OCT achieve a somewhat lower axial resolution.[8] The SS-Topcon device has a 100KHz A scan rate using a wavelength of 1050nm. A 3x3 volume scan was performed and each B-scans position was automatically scanned four times.[8] The Heidelberg prototype uses a full-spectrum amplitude decorrelation algorithm (FS-ADA) to evaluate motion contrast, which allows the evaluation of blood flow without sacrificing depth resolution. [9] The prototype acquired an A-scan rate of 85kHz with an axial resolution of 7 μ m and a lateral resolution of 14 μ m.[9] The available volume scan pattern was 4.3x 1.5 mm with 11 μ m between each B-scan. The ART frame was set at 13 frames per scan. Truetrack was employed to control for eye movements and minimize motion artefacts.

The study adhered to the tenets of the Declaration of Helsinki and was approved by the local ethics committee at Inselspital. A waiver of informed consent was granted due to the use of anonymized data (KEK number 2016–02100). The superficial (SCP) as well as deep vascular plexus (DCP) were segmented using the inbuilt software on each device. The location of the segmentation line of the SCP and the DCP of each module can be found in S1 Table in S1 Table and S2 Table file. Representative images for each module for the SCP and DCP are shown in Figs 1 and 2. The scans were checked for segmentation errors and subsequently the four separate, device-software generated en-face OCT-A images of the SCP and DCP from each volunteer were exported, analyzed and scored by 3 independent, experienced retinal imaging experts according to a pre-specified grading protocol.

Grading protocol

The grading protocol included qualitative as well as quantitative parameters: Motion artefacts (1 = no artefacts, 0 = some artefacts, -1 = severe motion artefacts), image artefacts (1 = no artefacts, 0 = some artefacts, -1 = severe image artefacts), the distinguishability of the foveal avascular zone (FAZ) (1 = FAZ border good distinguishable, 0 = middle, -1 = FAZ border barely/not distinguishable), and the vessel continuity and discriminability of large and small vessels (1 = vessel continuity clearly preserved, 0 = vessel continuity partly preserved, -1 = vessel continuity not preserved) were assessed. The number of clearly identifiable bifurcations identifiable on the superficial en-face image was counted on each scan. Therefore a main, large vessel

branch at 12 o'clock was chosen and the number of identifiable, subsequent bifurcations towards the terminal capillary end were counted on the respective branch. Additionally, the superficial layer retinal vessel density was assessed using the publicly available software Angio-tool. [10] Vessel density was defined as the area occupied by vessel lumens following binary reconstruction of images. [11].

Motion artifacts (or displacement artefacts) were considered as present when there were characteristic white-line artefacts present on OCT-A image with corresponding discontinuity at the en-face image in the B-scan direction or lateral displacement of parts of the image or doubling of retinal vessels. Axial motion artefacts resulting from breathing, tremor or pulsations as well as transverse motion artefacts caused by loss/change and saccades were evaluated. [12,13] Blink lines, identifiable as continuous dark lines of varying width visible on each singular slab were also included in this category. [14] Image artefacts were defined as any anomaly in the visual representation of information of the SCP and DCP slabs derived from the scanned object aside from motion artefacts. [12,14] These artefacts included segmentation artefacts, leading to deviation of the slab and projection artifacts, which were assumed when there were vessels seen clearly at deeper location than they actually inhabit. [12,15,16] It further included "negative projections" derived from superincumbent vessels after removing the projection flow signal using available device inbuilt software.

The parameter "some artefact" was defined as the presence of at least one artefact of respective category. The presence of 5 or more artefacts of respective category was considered as severe artefacts. Artefacts making reasonable evaluation of the microvasculature impossible were considered as severe artefacts as well. This included also broad artefacts as large as 5% of the image width or length in horizontal and vertical directions, respectively. Some representative examples of artifacts can be found in Figs 3 and 4.

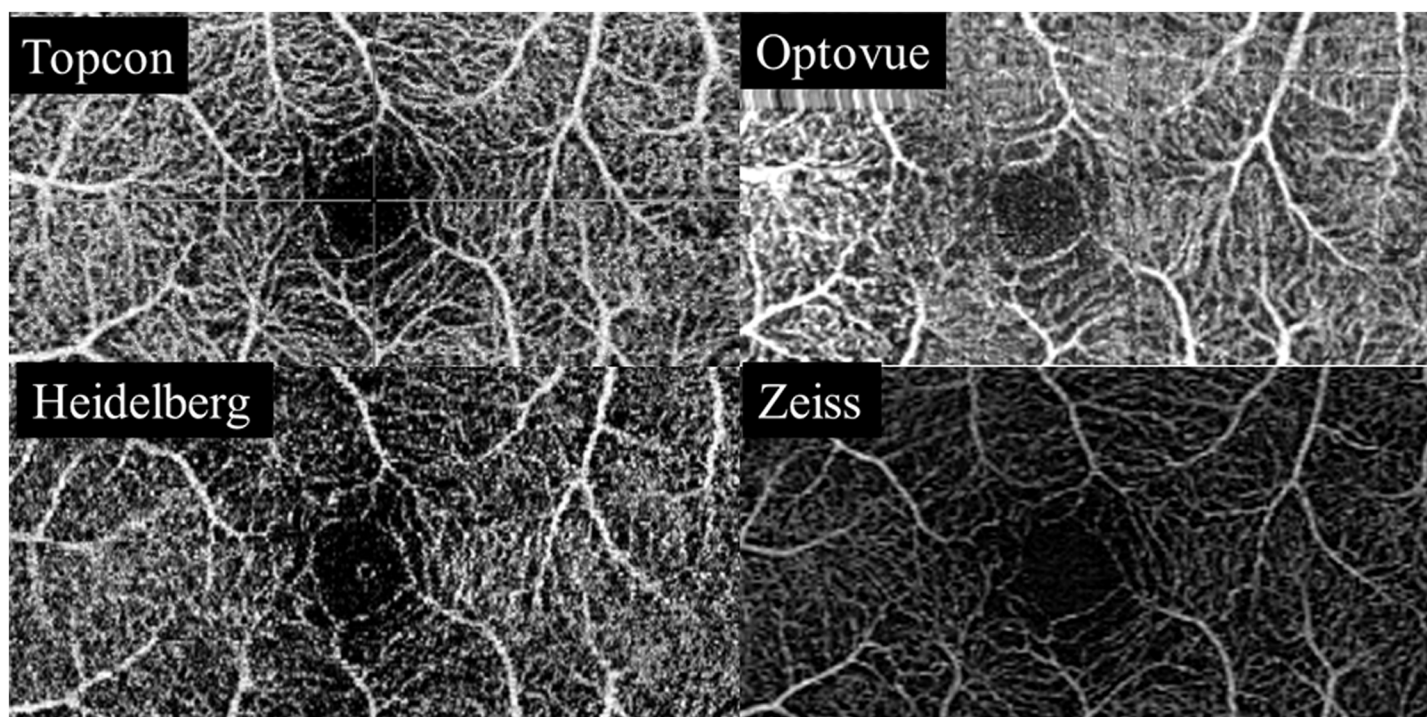


Fig 1. Superficial capillary plexus. Representative en face scans of the superficial capillary plexus (SCP) using the Swept source OCT Angio Topcon DRI OCT Triton (Top left), the Angiovue Optovue RTVue XR Avanti, (Top right), the Prototype of Spectralis OCT2 module with full spectrum decorrelation algorithm, Heidelberg Engineering (Bottom left) and the Zeiss AngioPlex Cirrus 5000 HD-OCT (Bottom right).

<https://doi.org/10.1371/journal.pone.0177059.g001>

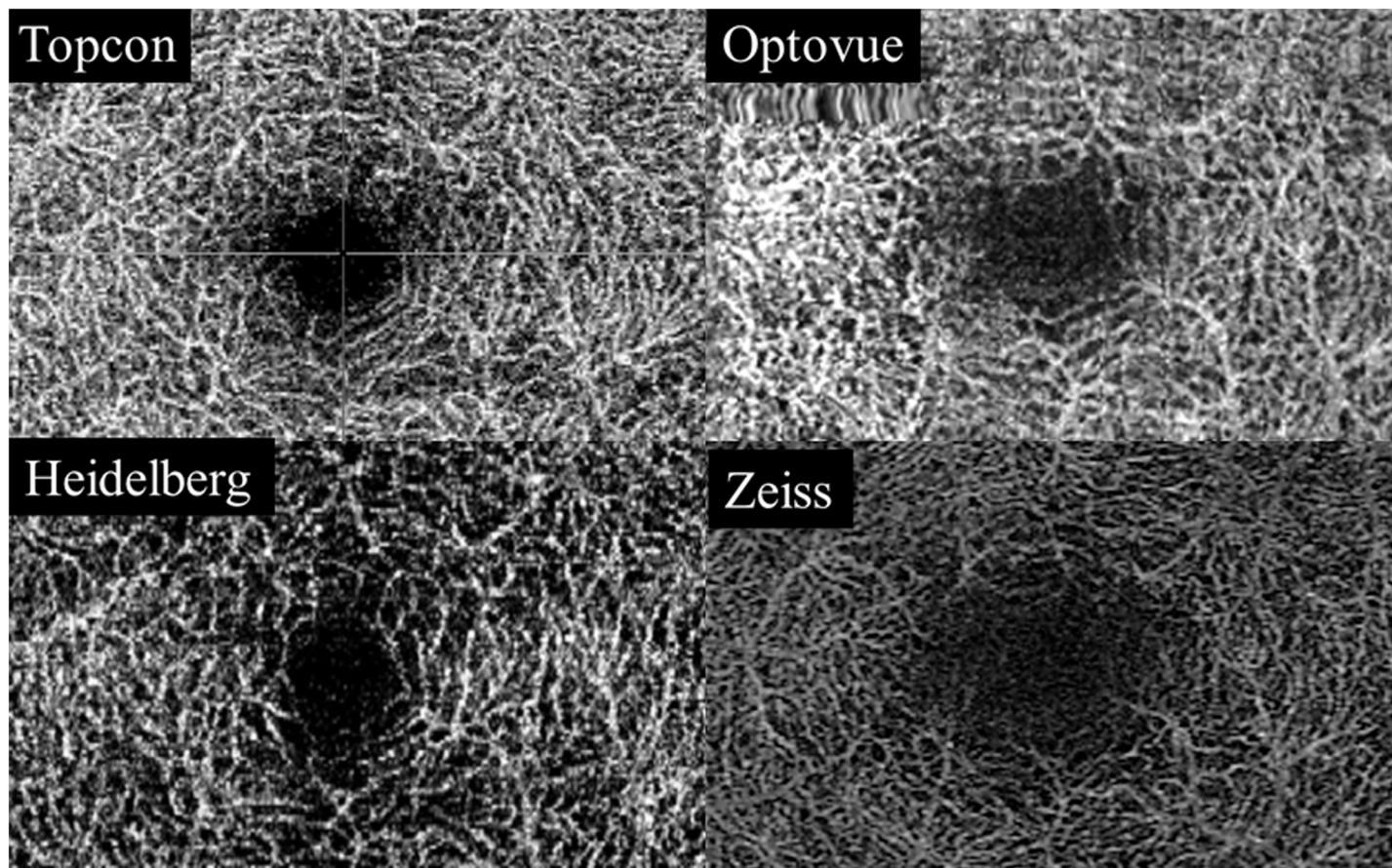


Fig 2. Deep capillary plexus. Representative en face scans of the deep capillary plexus (SCP) using the Swept source OCT Angio Topcon DRI OCT Triton (Top left), the Angiovue Optovue RTVue XR Avanti, (Top right), the Prototype of Spectralis OCT2 module with full spectrum decorrelation algorithm, Heidelberg Engineering (Bottom left) and the Zeiss AngioPlex Cirrus 5000 HD-OCT (Bottom right).

<https://doi.org/10.1371/journal.pone.0177059.g002>

After concordance analyses and the evaluation of the inter-grader reliability, a consensus grading was performed. The consensus dataset was then used to compare, score and rank the four modules.

Statistical analyses

Statistical analyses were performed with SPSS (IBM, SPSS statistics, Version 21, SPSS Inc., Chicago) software and R (www.r-project.org). Fleiss kappa coefficient was employed to quantify inter-grader reliability. After the evaluation of the interrater reliability, a consensus grading was performed and a final score ranging from -1/0/+1 was given for each feature for each image. These scores were then summed up and normalized with a maximal and minimal scoring ranging from -100 to +100. Based on the scores the devices were ranked for each evaluated parameter. Differences in ordinal variables such as motion artefacts, FAZ and vessel continuity were analyzed using Chi-squared test and numeric data such as number of counted bifurcations as well as vessel density were evaluated using ANOVA. An overall ranking of the modules was performed using exact binominal testing. Hereby it was evaluated in how many cases the module was better than the median with the set of evaluated parameters serving as the test sample for all features.

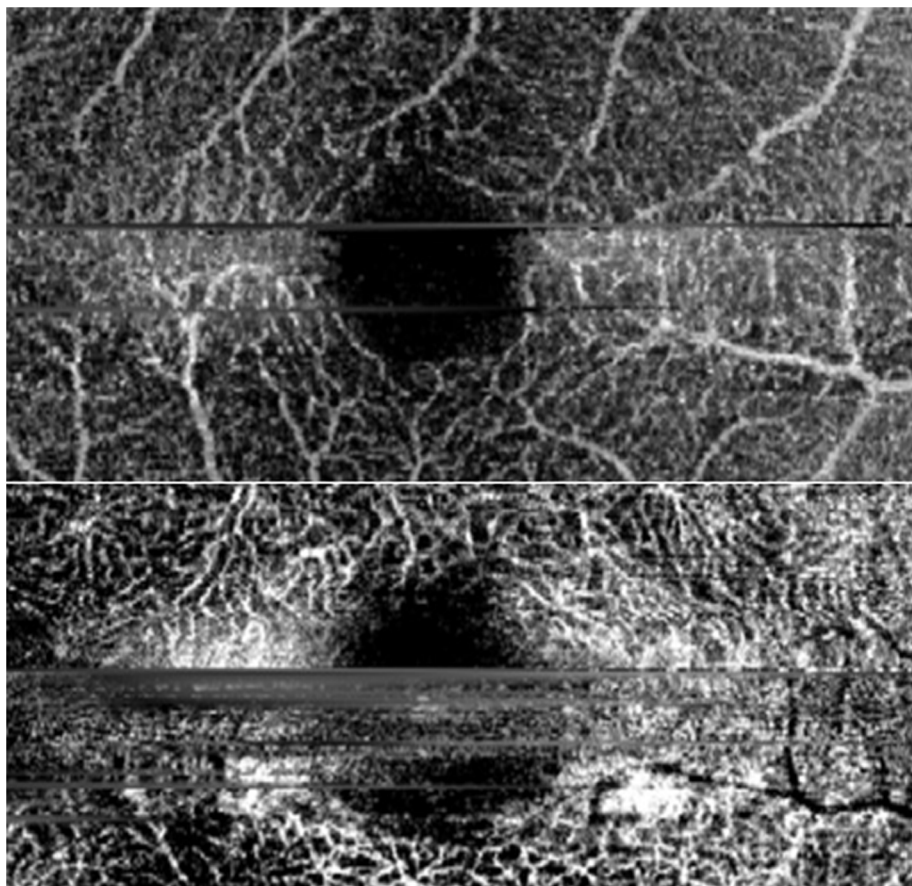


Fig 3. Severe motion artefacts. Representative examples of severe motion artefacts.

<https://doi.org/10.1371/journal.pone.0177059.g003>

P-values <0.05 were considered statistically significant and Bonferroni Holmes correction was used to adjust for multiple testing. Values are given as mean \pm SD.

Results

This study included 19 eyes of 19 healthy volunteers (mean age 35.3 ± 8.2 years). Inter-grader reliability in respect to each OCT-A module can be found in [Table 1](#).

OCT-A module differences

There was no difference in the overall vessel density among the evaluated modules using Angiotool (Zeiss $48.7\pm 4\%$, Optovue $47.9\pm 3\%$, Topcon $48.3\pm 2\%$, Heidelberg $46.5\pm 4\%$, $p = 0.2$). However although there was no difference in vessel density, the correlation coefficients were rather weak for respective parameter among the modules (Spearman correlation coefficient ranging from $r = -0.16$ – 0.35 , details see [S2 Table](#) in S1 and S2 Table file) No significant difference among the devices in terms of motion artefacts were detected ([Table 2](#)). However, for image artefacts of the SCP the Zeiss and the Topcon modules were superior compared to the other two devices ([Table 2](#)). The FAZ border of the SCP slabs were best appreciable on the Zeiss images, followed by the Optovue device, whereas the FAZ of the DCP was best discernable on the Optovue device followed by the Heidelberg module. The illustration of the continuity of the vessels was found to be superior on the Zeiss module ([Table 2](#)). The ranking of the

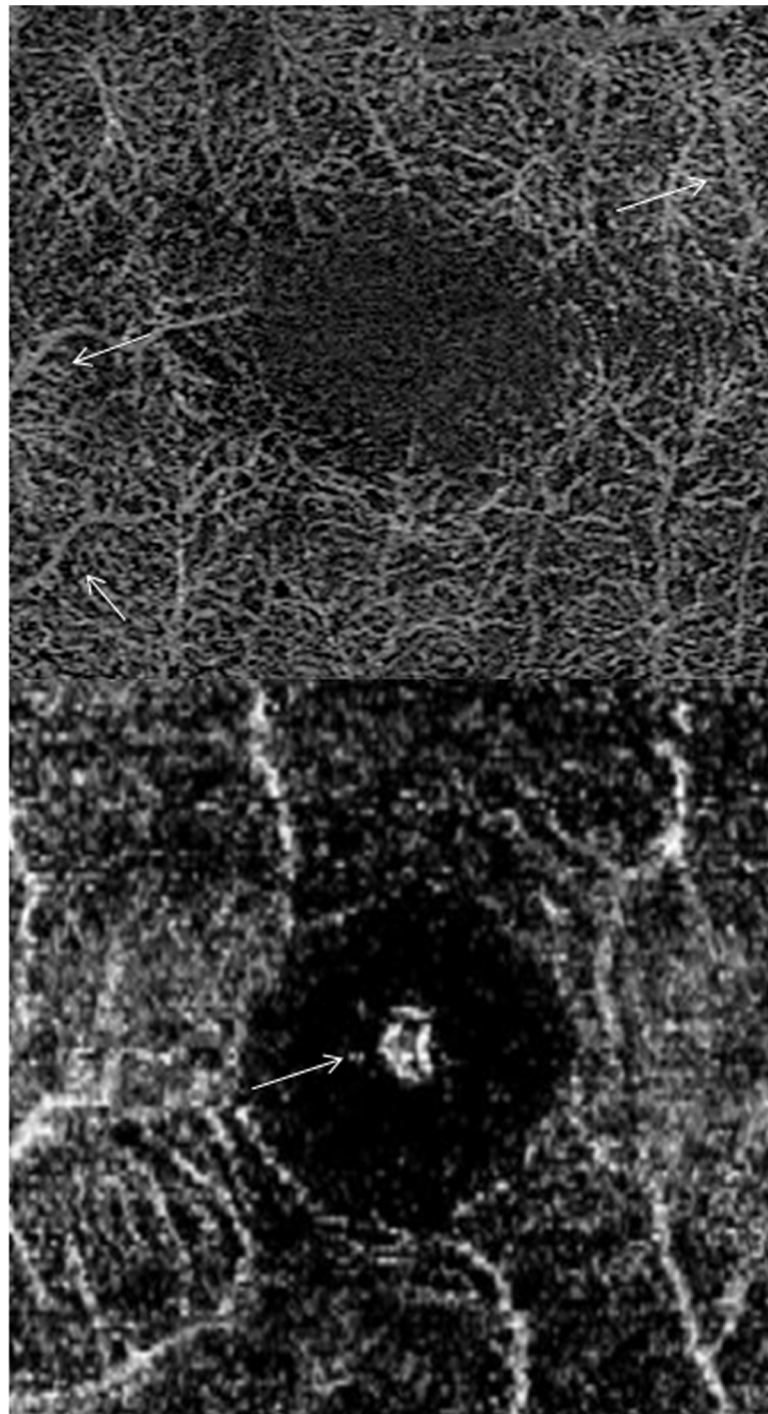


Fig 4. Projection artefacts and further examples of evaluated parameters. Top: Representative example of some projection artefacts (white arrows). Bottom: Example of “vessel continuity not preserved” of the small and large vessels. Central specular dots (white arrow) a form of image artefact, which can be seen in healthy eyes is shown.

<https://doi.org/10.1371/journal.pone.0177059.g004>

modules according to each individual evaluated feature can be found in [Table 2](#). The underlying normalized scores can also be found in [Table 2](#).

Table 1. Inter-grader reliability in respect to individual devices.

	Optovue		Topcon		Zeiss		Heidelberg	
	kappa	p-value	kappa	p-value	kappa	p-value	kappa	p-value
Motion artefacts SCP	0.818	<0.0001	0.566	<0.0001	0.527	<0.0001	0.31	≤0.002
Motion artefacts DCP	0.727	<0.0001	0.6	<0.0001	0.686	<0.0001	0.425	<0.0001
Image artefacts SCP	0.662	<0.0001	0.1	0.47	0.912	<0.0001	0.165	0.087
Image artefacts DCP	0.103	0.259	0.442	<0.0001	0.524	<0.0001	0.224	0.022
FAZ superior	0.515	<0.0001	0.396	<0.0001	0.782	<0.0001	0.31	0.00122
FAZ deep	0.535	<0.0001	0.107	0.367	0.493	<0.0001	0.357	0.002
Large vessel cont SCP	0.718	<0.0001	0.412	0.0002	0.912	<0.0001	0.285	0.007
Small vessel cont SCP	0.685	<0.0001	0.644	<0.0001	0.592	<0.0001	0.458	<0.0001
Small vessel cont DCP	0.613	<0.0001	0.245	0.0132	0.546	<0.0001	0.605	<0.0001
N of bifurcation	0.148	0.0453	0.295	0.0003	0.541	<0.0001	0.345	0.0006
Overall reliability SCP	0.522	<0.0001	0.296	0.0087	0.703	<0.0001	0.208	0.116
Overall reliability DCP	0.605	<0.0001	0.584	<0.0001	0.743	<0.0001	0.268	0.0133

Kappa values given in dark grey indicate a strong agreement (kappa value ranging from 0.9–0.7). Kappa values indicating a moderate agreement (= kappa value ranging from 0.5–0.7) are given in light grey. Kappa values of weak agreement (= kappa value ranging from 0.3–0.5) are shown in italic print style and minimal agreements (kappa value ranging from 0.1–0.3) are given in normal print style. Overall reliability indicate the intergrader reliability of all evaluated features in the SCP and the DCP together. Kappa = Fleiss kappa, SCP = superficial capillary plexus, DCP = deep capillary plexus, FAZ = foveal avascular zone, cont. = continuity, N = number.

<https://doi.org/10.1371/journal.pone.0177059.t001>

Significant different numbers of bifurcations were discernible on each module (Zeiss 2 ± 0.9 bifurcations, Optovue 2.5 ± 1.2 , Topcon 1.3 ± 0.7 and Heidelberg 0.5 ± 0.6 , $p \leq 0.001$, Table 2). In the overall ranking, the Zeiss module was superior and in 90% better than the median (Bonferroni corrected p-value = 0.04, Fig 5). The Optovue was found to be better than the median in

Table 2. Ranking and underlying normalized scores of each module for each evaluated variable of the consensus dataset.

	Ranking (Scores)				p-value
	Optovue	Topcon	Zeiss	Heidelberg	
Motion artefacts SCP	2 (42)	3 (26)	1 (74)	4 (11)	0.135
Motion artefacts DCP	2 (42)	3 (21)	1 (63)	4 (15)	0.076
Image artefacts SCP	3 (58)	2 (74)	1 (74)	4 (-21)	≤0.001
Image artefacts DCP	4 (12)	1 (30)	2 (19)	3 (14)	0.32
FAZ SCP	2 (47)	3 (36)	1 (84)	4 (5)	≤0.001
FAZ DCP	1 (31)	3 (-37)	4 (-58)	2 (-21)	0.002
Large vessel cont SCP	3 (47)	2 (63)	1 (74)	4 (26)	0.027
Small vessel cont SCP	2 (53)	3 (16)	1 (74)	4 (-58)	≤0.001
Small vessel cont DCP	3 (36.7)	2 (37.1)	1 (42)	4 (-11)	0.071
N of bifurcation	1 (2.5±1.2)	3 (1.3±0.7)	2 (2±0.9)	4 (0.5±0.6)	≤0.001

Each image was graded by the readers and was given a score ranging from -1/0/+1. (e.g. Motion artefacts (1 = no artefacts, 0 = some artefacts, -1 = severe motion artefacts). Thereafter a consensus grading was performed and a final score was given for each feature for each image. These scores were then summed up and normalized. The maximal and minimal possible scoring after normalization was -100 to +100 (see numbers in brackets). Based on these scores each device was ranked (bold numbers). Rank of 1 describes highest scores and best quality (in grey). The variable “n of Bifurcation” shows mean and standard deviation of counted bifurcations (see numbers in brackets). Therefore the main, large vessel branch at 12 o'clock was chosen and the number of identifiable, subsequent bifurcations towards the terminal capillary end were counted on the respective branch and respective variable was then ranked based on the number of identifiable bifurcations (bold numbers). Differences among the modules (p-values) were calculated using Chi-squared testing.

<https://doi.org/10.1371/journal.pone.0177059.t002>

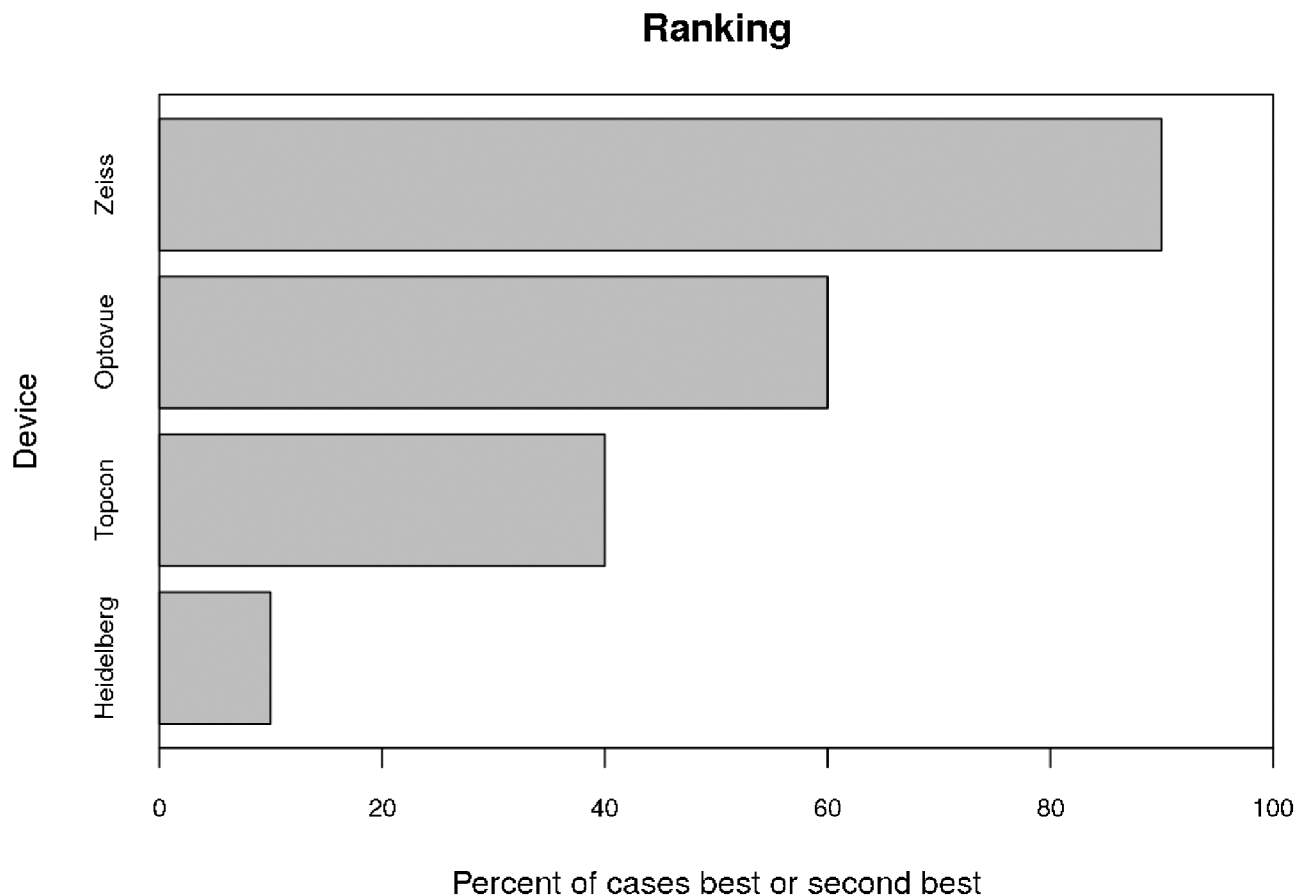


Fig 5. OCT-A module ranking.

<https://doi.org/10.1371/journal.pone.0177059.g005>

60%, the Topcon in 40% and the Heidelberg module in 10%, however these evaluations missed statistical significance (Fig 5).

Overall ranking of each OCT-A module using exact binominal testing based on the scores of each evaluated feature. In 90% the Zeiss Angioplex module was better than the median (corrected p-value 0.04), The Optovue was in 60%, the Topcon in 30% and the Heidelberg in 10% better than the median. However, differences are not statistically significant.

Discussion

We here present the first comparison of 4 OCT-A modules. Each of the four modules had particular strengths, which differentiated it from their competitors. Overall, however, the Zeiss module seemed to be the most reliable, accurate and precise device in terms of our evaluated variables, followed by the Optovue, the Topcon and the Heidelberg module. Each module employed different technology to quantify the motion contrast and each module had different approaches to minimize motion artefacts and achieve optimal image quality with high resolution. [6–8,17].

Inter-grader agreement also differed in respect to evaluated feature and device. It seems noteworthy that the reliability of the grading was found to be higher in the Zeiss and Optovue module compared to the Topcon and Heidelberg module. This may be explained by the fact that the grader reliability seems associated with the quality of the images and evaluated features. Images of higher quality depicting the evaluated feature more accurately will result in a

stronger inter-grader agreement than poor images. Taking into account that overall the Zeiss module was ranked best and the Optovue device second best, the differences in the inter-grader reliability seem reasonable.

Motion artefacts in the SCP and the DCP were less prominent with the Zeiss Angioplex module, followed by the OCT-A modules of Optovue, Topcon and Heidelberg. However, these findings were not statistically significant. A previous study compared the AngioVue (Optovue) with the Angioplex (Zeiss) and found that the Angioplex required shorter acquisition time and showed a lower number of motion artefacts when compared to the Angiovue. [13] Further the number of low signal strength images and the images impossible to analyze were lower in the Angioplex module compared to the Angiovue module.[13].

In order to prevent motion artefacts manufacturers use different approaches. Zeiss (FastTrack), Heidelberg (TrueTrack) and Topcon (SMARTTRACK) implemented a retinal eye tracking, while the here evaluated Optovue device, used a software based method in which a retinal area is repeatedly scanned horizontally and vertically. This software based approach may account for longer acquisition times compared to the eye tracking as shown in the previous study by De Vities et al. [13] In our study, however it did not seem to impact the severity of motion artefacts as the Optovue module was ranged second best for this evaluated feature. Further, the now available Angiovue modules provide real time eye tracking as well. It has been claimed that the SSADA algorithm mainly accounts for axial artefacts and therefore transverse artefacts may still remain an issue. This assumption was also not confirmed in our analyses. Although not statistically significant, the Zeiss module scored highest in terms of the absence of motion artefacts. This may be explained by the fact that beside the Fasttrack technology for continuous eye tracking, this module also samples the retina 15 times per second to minimize motion artefacts. Only areas which may be affected by motion artefacts are rescanned which reduces the acquisition time and thereby again motion artefacts. [6].

To sum up, the presence and severity of motion artefacts did not significantly differ among the evaluated modules and manufacturers are continuously working on better solutions to delimitate respective problems on OCT-A. For instance, some manufacturers now offer motion correction technology in order to overcome respective artefacts.

The category of imaging artefacts in our study comprised the presence of segmentation and projection artifacts. Projection artefacts arise from light, which is not directly reflected by the moving blood but passes through and illuminate features posterior to the vessel. [12] This implicates that mainly the slabs of the DCP were affected by respective artefacts in this study, while segmentation artefacts can be found in the SCP and DCP. Projection artefacts occur in all quantifying motion contrast methods irrespective whether speckle-or intensity decorrelation or phase variance is used. [18] Many manufacturers now offer software implementation which remove respective artifacts.[19] However, after removal of the projection flow signal using device inbuilt software, a “negative projection” visible as dark shadow of the same vascular pattern remains. Zeiss, Topcon as well as Optovue provided such inbuilt image processing for the removal of projection artefacts, whereas the Heidelberg prototype did not have respective post acquisition processing available at the timepoint of our evaluation. In terms of image artefacts the Zeiss and Topcon was superior to the Optovue and Heidelberg module, although just the alterations found in the SCP were statistical significant. One reason why Topcon was superior to Optovue and Heidelberg may be that SS- OCT at a wavelength of 1050 nm was used. Longer wavelengths are less susceptible to light scattering, which may decrease the presence of projection artefacts in the images of the DCP. The Angioplex uses the so-called OCT-microangiography complex (OMAG), which identifies changes in the phase and intensity information of the OCT scans. [7] The reason why the Angioplex was superior compared to the Optovue and the Heidelberg is not explainable by a longer wavelength. The use of OMAG

together with the above-described advantageous implementations of the Zeiss module may account for this fact.

A recent study suggested that the usage of the FS-ADA algorithm, which is implemented in the Heidelberg module, would cause less projection artefacts compared to the SSADA algorithm.[20] The SSADA algorithm creates an isotropic voxel by degrading the axial resolution until the axial and the transverse dimension is equal, which can cause projection artefacts. In contrast, the FS-ADA algorithm detects flow from structural OCT images without impairing the axial resolution.[20] On the DCP en-face images of our analyses the Heidelberg module using FS-ADA showed indeed numerically less severe image artifacts than the Optovue, however on the SCP the Heidelberg module was inferior. This may be in line with the previous assumptions because projection artefacts are more likely to be found on the deeper en-face slabs. [12,20] Of course the inbuilt software for projection artefact removal in the Zeiss, Optovue and Topcon modules have probably also led to superior results, although in many cases “negative projections” were seen in respective cases.

But these explanations only account for the artefacts due to projections in the DCP. For the SCP mainly segmentation errors attribute to the artefacts in the category. Previous studies have shown that segmentation artefacts are rather common. [16,21] For layer segmentation different approaches are applied in each module including preprocessing steps for OCT denoising and methods such as pattern recognition, pixel classification of retinal layers, graph based multi-surface segmentation, global segmentation algorithms including active contours and Markov random fields, artificial intelligence approaches based on multiresolution hierarchical support vector machine or fuzzy C-means clustering techniques and 3D graph based methods, which lead to more or less accurate retinal layer segmentation.[22] An extensive discussion of the pros and cons of the different segmentation approaches used definitely exceeds the scope of this paper.

Interestingly the distinguishability of the FAZ borders differed between the SCP and the DCP. While the FAZ borders of the SCP was best visualized by the Zeiss Angioplex followed by the Optovue, Topcon and Heidelberg, the borders of the DCP were harder to distinguish and were best identifiable on the Optovue followed by the Heidelberg, Topcon and Zeiss module. The SSADA algorithm employed on the Optovue uses a four-fold spectrum split, which improves the signal to noise ratio, which may account for the high scoring on the distinguishability of the FAZ border on the SCP and the DCP. [7,18] It was previously shown that the SSADA algorithm provides a clean and continuous microvascular network and barely noise inside the FAZ. [7] This observation can be confirmed by our observation as the Optovue was superior in distinguishing the FAZ border in the DCP and was the second best module in differentiating the FAZ in the SCP. The ART frame on the Heidelberg device, which was set at 13 frames per B-scan, may have enhanced the signal to noise ratio and may have therefore improved the visibility of the FAZ of the DCP, which was in general harder to distinguish than the FAZ of the SCP (see scoring in Table 2). The rather inferior performance of the Topcon module in terms of respective parameter may be partly explained by the usage of a 1050nm wavelength, because spectral domain OCT using a wavelength around 800nm produces high quality angiograms with less axial scans needed and more transverse points in less time. This may be due to a lower decorrelation noise needing only two consecutive scans instead of eight for one position. [7,18] Another reason may be the different segmentation boundaries of each device. [23,24] Spaide et al. found that the different segmentation boundaries and their default segmentation result in different sizes and shapes of the FAZ. [24] These differences mainly originate from the inner retinal layers, which become thinner and terminate as they reach the central fovea.[24] Another paper suggested that due to the so far incorrect anatomical segmentation algorithms and the great inter-individual disparity, the best approach for the

visualization of the FAZ would be the usage of the whole retina slabs. [23] Beside the above mentioned factors, the different approaches to quantify motion contrast in each module, the eye tracking and the methods to increase resolution and signal to noise ratio, other factors such as predefined contrast settings may have impacted the discrimination of the FAZ borders on each device.[25].

Vessel continuity on the SCP and DCP were best preserved and discriminable on the Zeiss module, followed by the Topcon, Optovue and Heidelberg. Of course, the severity of motion and imaging artefacts had significant impact on this evaluated parameter as high resolution and the absence of artefacts is key for the continuity and discriminability of vessels. Therefore already above mentioned parameters may account for the superior presentability of the vessels on the Zeiss module. A previous investigation indicated that the vascular network may be better visualized using OCTARA (Topcon) than the SSADA module. [8] It is further known that the SSADA algorithm mainly accounts for axial artefacts, whereas transverse artefacts, which may also cause discontinuity of vessels, remain a problem. These findings may be in line with our observations as the Topcon module was second best for the discernibility of the large vessels on the SCP and the small vessels on the DCP. However, small vessels on the SCP were better identifiable using Optovue compared to Topcon.

The vessel density of all 4 devices was evaluated using the publicly available Angiotool and was comparable among all 4 modules. The measured vessel density of the SCP in our cohort of all devices was also comparable to the normative vessel density previously reported by Coscas et al., which was evaluated with the inbuilt Angiovue software and found to be around 52.58 \pm 3.22%. [26] However although there was no difference in means, the correlation among the devices was rather weak, implicating that artefacts and differences in terms of the FAZ size and contour significantly impact the evaluation of the vessel density. This important fact should be considered when evaluating the vessel density and special attention should be drawn to the quality of the evaluated scans.

This study has definite limitations. First, a prototype was tested against three commercially available modules and it remains to be shown whether the final version of the Heidelberg prototype improves their performance compared to the other so far available modules. Second, the scanned area on the Heidelberg device (4.3x 1.5 mm) differed from the area scanned with the remaining devices (3x3mm). This may have impacted resolution and acquisition time and limits the comparability to the other devices. Usually image quality of healthy subjects outperforms quality of images of diseased eyes. Thus, the evaluation and comparison of respective modules in diseased eyes are warranted as well. Beside the 4 evaluated modules there are other OCT-A systems currently under development which were not tested in this study, including the OCT RS-3000Advance of Nidek and the SOCT Copernicus REVO and REVO NX of Optopol.

In conclusion, each device uses different approaches to offer optimal high-resolution images of the vascular network with the minimum possible artefacts. In our study, the approach of Zeiss was in most of the evaluated features superior and better than the median, but all modules had their individual strengths and weaknesses. The current study reflects just the current stage of development, but the OCT-A imaging is still in its early beginnings and is rapidly developing and improving with great strides forward.

Supporting information

S1 Table. Location of boundaries of each evaluated device outlining the superficial capillary plexus (SCP) and the deep capillary plexus (DCP).
(DOCX)

S2 Table. Spearman correlation coefficients among the modules in terms of vessel density. (DOCX)

Acknowledgments

The authors would like to thank Rene Rückert, MD, MBA for the critical discussion and review of the manuscript.

Author Contributions

Conceptualization: MRM SW MSZ.

Data curation: MRM HGZ LB WH AE.

Formal analysis: WH MRM.

Funding acquisition: SW MSZ.

Investigation: MRM HGZ LB.

Methodology: MRM SW MSZ.

Project administration: MRM.

Resources: SW MSZ AE.

Software: WH AE.

Supervision: MRM MSZ SW.

Validation: WH MRM HGZ LB.

Visualization: MRM.

Writing – original draft: MRM.

Writing – review & editing: MRM HGZ LB WH AE SW MSZ.

References

1. Lipson BK, Yannuzzi LA Complications of intravenous fluorescein injections. *Int Ophthalmol Clin*. 1989; 29:200–205. PMID: [2526795](#)
2. Park JJ, Soetikno BT, Fawzi AA Characterization of the Middle Capillary Plexus Using Optical Coherence Tomography Angiography in Healthy and Diabetic Eyes. *Retina*. 2016.
3. Zhang A, Zhang Q, Chen CL, Wang RK Methods and algorithms for optical coherence tomography-based angiography: a review and comparison. *J Biomed Opt*. 2015; 20:100901. <https://doi.org/10.1117/1.JBO.20.10.100901> PMID: [26473588](#)
4. Gao SS, Jia Y, Zhang M, Su JP, Liu G, Hwang TS, et al. Optical Coherence Tomography Angiography. *Invest Ophthalmol Vis Sci*. 2016; 57:OCT27–36. <https://doi.org/10.1167/iovs.15-19043> PMID: [27409483](#)
5. Ruminski D, Sikorski BL, Bukowska D, Szkulmowski M, Krawiec K, Malukiewicz G, et al. OCT angiography by absolute intensity difference applied to normal and diseased human retinas. *Biomed Opt Express*. 2015; 6:2738–2754. <https://doi.org/10.1364/BOE.6.002738> PMID: [26309740](#)
6. Rosenfeld PJ, Durbin MK, Roisman L, Zheng F, Miller A, Robbins G, et al. ZEISS Angioplex Spectral Domain Optical Coherence Tomography Angiography: Technical Aspects. *Dev Ophthalmol*. 2016; 56:18–29. <https://doi.org/10.1159/000442773> PMID: [27023249](#)
7. Huang D, Jia Y, Gao SS, Lumbroso B, Rispoli M Optical Coherence Tomography Angiography Using the Optovue Device. *Dev Ophthalmol*. 2016; 56:6–12. <https://doi.org/10.1159/000442770> PMID: [27022989](#)

8. Stanga PE, Tsamis E, Papayannis A, Stringa F, Cole T, Jalil A Swept-Source Optical Coherence Tomography Angio (Topcon Corp, Japan): Technology Review. *Dev Ophthalmol*. 2016; 56:13–17. <https://doi.org/10.1159/000442771> PMID: 27023108
9. Coscas G, Lupidi M, Coscas F Heidelberg Spectralis Optical Coherence Tomography Angiography: Technical Aspects. *Dev Ophthalmol*. 2016; 56:1–5. <https://doi.org/10.1159/000442768> PMID: 27022921
10. Zudaire E, Gambardella L, Kurcz C, Vermeren S A computational tool for quantitative analysis of vascular networks. *PLoS One*. 2011; 6:e27385. <https://doi.org/10.1371/journal.pone.0027385> PMID: 22110636
11. Coscas G, Coscas F, Zucchiatti I, Glacet-Bernard A, Soubrane G, Souied E SD-OCT pattern of retinal venous occlusion with cystoid macular edema treated with Ozurdex(R). *Eur J Ophthalmol*. 2011; 21:631–636. <https://doi.org/10.5301/EJO.2011.7428> PMID: 21500185
12. Spaide RF, Fujimoto JG, Waheed NK Image Artifacts in Optical Coherence Tomography Angiography. *Retina*. 2015; 35:2163–2180. <https://doi.org/10.1097/IAE.0000000000000765> PMID: 26428607
13. De Vitis LA, Benatti L, Tomasso L, Baldin G, Carnevali A, Querques L, et al. Comparison of the Performance of Two Different Spectral-Domain Optical Coherence Tomography Angiography Devices in Clinical Practice. *Ophthalmic Res*. 2016.
14. Say EA, Ferenczy S, Magrath GN, Samara WA, Khoo CT, Shields CL IMAGE QUALITY AND ARTIFACTS ON OPTICAL COHERENCE TOMOGRAPHY ANGIOGRAPHY: Comparison of Pathologic and Paired Fellow Eyes in 65 Patients With Unilateral Choroidal Melanoma Treated With Plaque Radiotherapy. *Retina*. 2016.
15. Ferrara D Image artifacts in optical coherence tomography angiography. *Clin Experiment Ophthalmol*. 2016; 44:367–368. <https://doi.org/10.1111/ceo.12781> PMID: 27381573
16. Ghasemi Falavarjani K, Al-Sheikh M, Akil H, Sadda SR Image artefacts in swept-source optical coherence tomography angiography. *Br J Ophthalmol*. 2016.
17. Lumbroso B, Rispoli M, Savastano MC, Jia Y, Tan O, Huang D Optical Coherence Tomography Angiography Study of Choroidal Neovascularization Early Response after Treatment. *Dev Ophthalmol*. 2016; 56:77–85. <https://doi.org/10.1159/000442782> PMID: 27022967
18. Jia Y, Tan O, Tokayer J, Potsaid B, Wang Y, Liu JJ, et al. Split-spectrum amplitude-decorrelation angiography with optical coherence tomography. *Opt Express*. 2012; 20:4710–4725. <https://doi.org/10.1364/OE.20.004710> PMID: 22418228
19. Zhang Q, Zhang A, Lee CS, Lee AY, Rezaei KA, Roisman L, et al. Projection Artifact Removal Improves Visualization and Quantitation of Macular Neovascularization Imaged by Optical Coherence Tomography Angiography *Ophthalmology Retina*. 2016; epub ahead of print.
20. Lupidi M, Coscas F, Cagini C, Fiore T, Spaccini E, Fruttini D, et al. Automated Quantitative Analysis of Retinal Microvasculature in Normal Eyes on Optical Coherence Tomography Angiography. *Am J Ophthalmol*. 2016; 169:9–23. <https://doi.org/10.1016/j.ajo.2016.06.008> PMID: 27296485
21. Alshareef RA, Dumpala S, Rapole S, Januwada M, Goud A, Peguda HK, et al. Prevalence and Distribution of Segmentation Errors in Macular Ganglion Cell Analysis of Healthy Eyes Using Cirrus HD-OCT. *PLoS One*. 2016; 11:e0155319. <https://doi.org/10.1371/journal.pone.0155319> PMID: 27191396
22. Kafieh R, Rabbani H, Kermani S A review of algorithms for segmentation of optical coherence tomography from retina. *J Med Signals Sens*. 2013; 3:45–60. PMID: 24083137
23. Campbell JP, Zhang M, Hwang TS, Bailey ST, Wilson DJ, Jia Y, et al. Detailed Vascular Anatomy of the Human Retina by Projection-Resolved Optical Coherence Tomography Angiography. *Sci Rep*. 2017; 7:42201. <https://doi.org/10.1038/srep42201> PMID: 28186181
24. Spaide RF, Curcio CA Evaluation of Segmentation of the Superficial and Deep Vascular Layers of the Retina by Optical Coherence Tomography Angiography Instruments in Normal Eyes. *JAMA Ophthalmol*. 2017; 135:259–262. <https://doi.org/10.1001/jamaophthalmol.2016.5327> PMID: 28097291
25. Palma CV, Amin R, Huf W, Schlanitz F, Eibenberger K, Jampol LM, et al. Spectral Domain-Optical Coherence Tomography Image Contrast and Background Color Settings Influence Identification of Retinal Structures. *Retina*. 2016; 36:1888–1896. <https://doi.org/10.1097/IAE.0000000000001060> PMID: 27219667
26. Coscas F, Sellam A, Glacet-Bernard A, Jung C, Goudot M, Miere A, et al. Normative Data for Vascular Density in Superficial and Deep Capillary Plexuses of Healthy Adults Assessed by Optical Coherence Tomography Angiography. *Invest Ophthalmol Vis Sci*. 2016; 57:OCT211–223. <https://doi.org/10.1167/iovs.15-18793> PMID: 27409475